

10 19

GSI Medical Documentation Pig

Uniform

Time of arrival at GSI:

10:53

Weight [kg]

39.1

Vital functions (arrival, in cave, immobilization, positioning, field #1, field #2)

| | | | | | | |
|------|-------|-------|-------|-------|--|--|
| TIME | 11:08 | 11:20 | 12:10 | 12:22 | | |
| HR | 47 | 34 | 66 | 66 | | |
| Temp | | | | | | |
| SpO2 | 96 | 88 | 100 | 99 | | |

Premedication:

Ketamin [mg] 100

Diazepam [mg] 15

Azaperon [mg] 160

+ 160 (11:01)

Preparation:

- ☒ Place IV in ear vein ☒ Intubate pig ID 7.0 ET
☐ Shave ECG areas ☒ Start 0.9% saline iv.
☒ Start propofol [ml/KG/min] 0.19

11:05 Propofol

11:20 80 mg Propofol

0.49 mg

Drawing of venous blood samples before immobilization

For immediate centrifugation @ 2800 rpm:

- ☒ 2x for Citrate ☒ 2x EDTA ☒ 2x Li-Hep

For centrifugation after 30 min @ 2800 rpm:

- ☒ 2x Serum

Transfer to cave

- ☒ Establish ECG monitoring: 3 ECG surface leads at shaved locations

Respiration:

- ☒ 14 l/min, IPPV, $V_t = 350$, PEEP 0 ☐ other: _____

Irradiation**Drawing of venous blood samples after irradiation**

For immediate centrifugation @ 2800 rpm:

- ☐ 2x for Citrate ☐ 2x EDTA ☐ 2x Li-Hep

For centrifugation after 30 min @ 2800 rpm:

- ☐ 2x Serum

Recovery:

Monitor SpO2 and HR until awake

Pig awake and healthy: ☐ Yes☐ No

Time: _____

Comments:

PIG: Uniform



Important Treatment Planning Parameters:


| | | CHECK |
|-------------------------------|-----------------------------------|-------|
| Dose: | 40 Gy | |
| Target: | AblPV | |
| Isocenter AblPVMargin5mm: | 256.0 217.0 71.0 | |
| BB Center: | 252.6 255.5 57.0 | |
| Translation vector: | 3.4 -38.5 -14.0 | |
| Final Table Position: | | |
| RST field 1: | RESCANNING_Uniform_90.rst | |
| Sam file produced? | | |
| Angle: | 90 | ✓ |
| Number of EOPs: | 19 | ✓ |
| Minimal energy: | 164.63 | |
| Maximal energy: | 237.05 | |
| Total number of rasterpoints: | 17824 | ✓ |
| Rasterpoints first plane: | 53 | ✓ |
| Rasterpoints last plane: | 2205 | ✓ |
| RST field 2: | RESCANNING_Uniform_270.rst | |
| Sam file produced? | | |
| Angle: | 270 | |
| Number of EOPs: | 22 | ✓ |
| Minimal energy: | 143.79 | |
| Maximal energy: | 232.20 | |
| Total number of rasterpoints: | 17526 | ✓ |
| Rasterpoints first plane: | 42 | ✓ |
| Rasterpoints last plane: | 1590 | ✓ |

12:33 -
12:49
16:49 min


13:25 -
13:43
17:43 min

(MPT: max 0.7 weight 0.9)

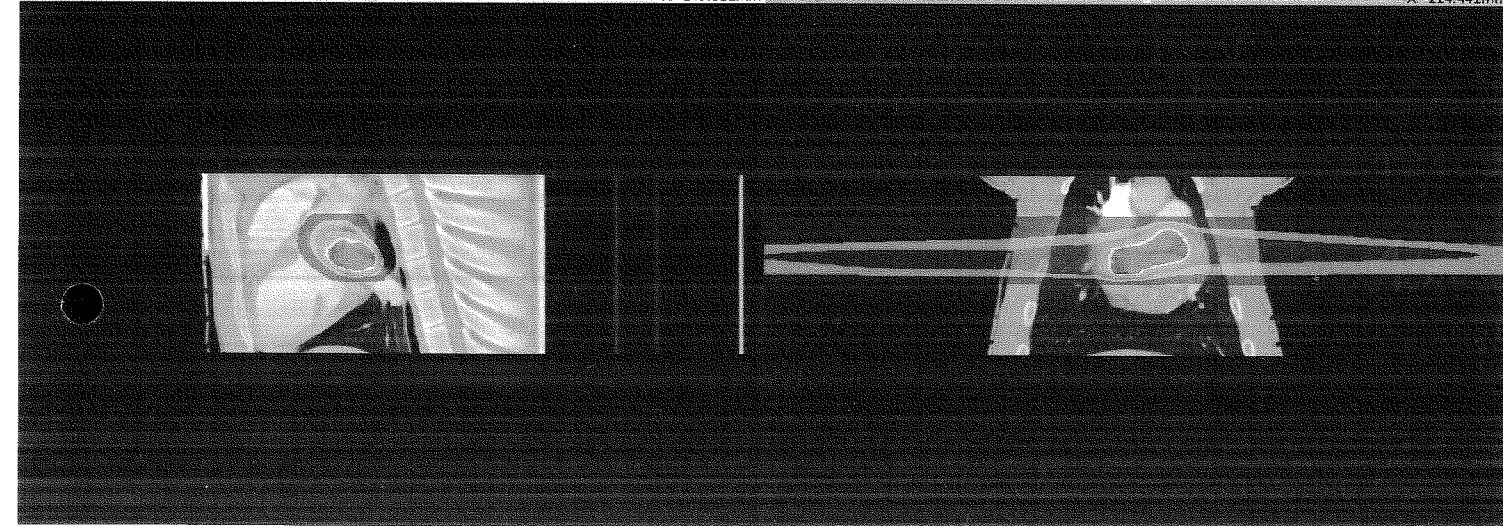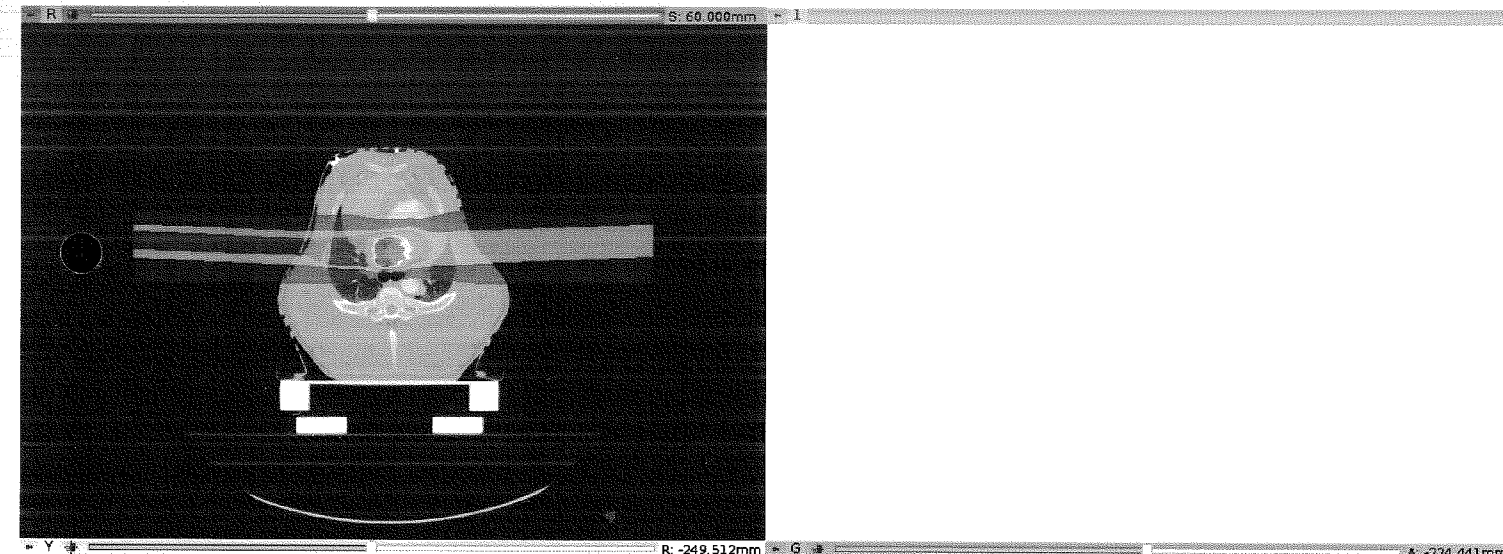

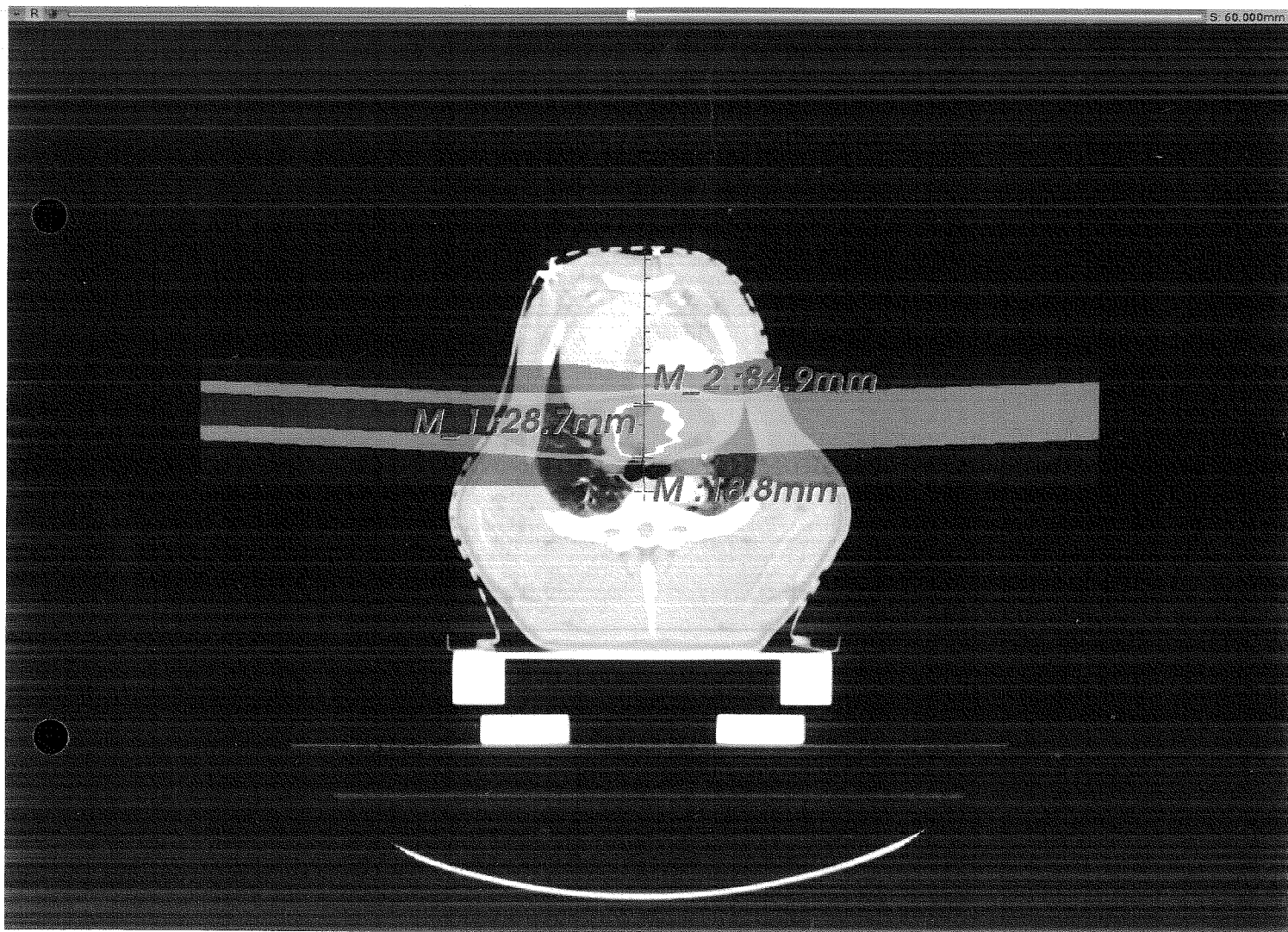

Uniform Abstände

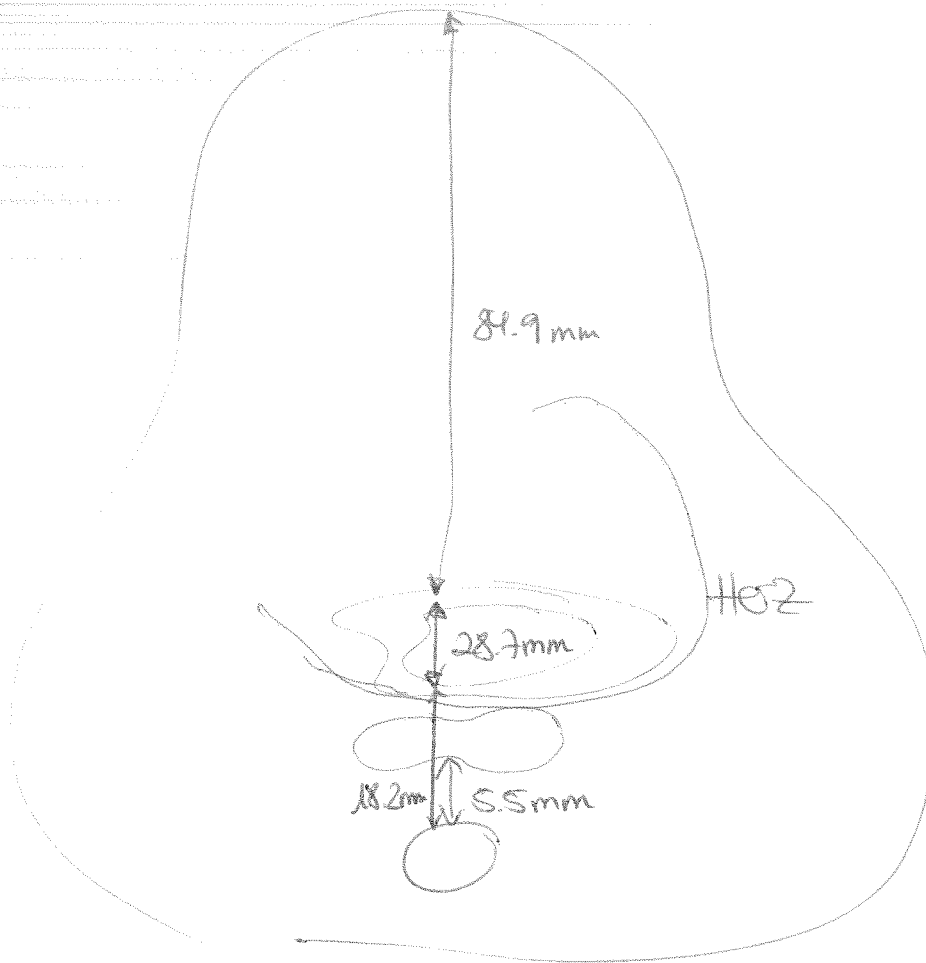

hinterranke Herz - Trachea Abstand

Positioning Documentation PigUniformImmobilization: PiggyFix ☐ A ☐ B Mask left: 4-8 right: 4-8☒ Match Tatoos, mask and laser, Beekleys in isocenter – check images on page #2!☐ patient couch zeroed ☐ Translation to target: X 3.4 Y -38.5 Z -14.0

Comments: _____

Respiration: ☒ 14 bpm, IPPV, $V_t = 350$, PEEP 0 ☐ other: _____Vertical X-ray #1 Time (HH:MM) 11:43 breathhold: _____ secsExposure: 77 kV 50 mA 1000 ms☒ Detector position ☒ Trigger cable set ☒ XIS external triggerX-ray filename Uniform-Vertical1 ☐ XIS – file closedLateral X-ray #1 Time (HH:MM) 11:40 breathhold: _____ secsExposure: 77 kV 50 mA 1000 ms☒ Detector position ☒ Trigger cable set ☒ XIS external triggerX-ray filename Uniform Lateral1 ☐ XIS – file closedVector: X 6.4 Y 7.1 Z 3.0 Position: X 9.8 Y -31.4 Z -10.4Match on: ☐ bony anatomy ☐ markers ☐ pacemaker ☒ bony + heart

Repeat x-rays on separate page! Number of x-rays _____

Comments: _____

Final vertical X-ray Time (HH:MM) _____ breathhold: _____ secs

Exposure: _____ kV _____ mA _____ ms

☐ Detector position ☐ Trigger cable set ☐ XIS external triggerX-ray filename _____ ☐ XIS – file closed

Final lateral X-ray Time (HH:MM) _____ breathhold: _____ secs

Exposure: _____ kV _____ mA _____ ms

☐ Detector position ☐ Trigger cable set ☐ XIS external triggerX-ray filename _____ ☐ XIS – file closed

Position adjustment: X _____ Y _____ Z _____

Final position: X 9.8 Y -31.4 Z -10.4☐ Pictures left/right☐ X-rays back-up

Additional X-rays Pig Uniform

Vertical X-ray #2

Time (HH:MM) 11:53

breathhold: _____ secs

Exposure: 77 kV 50 mA 1000 ms☒ Detector position ☒ Trigger cable set ☒ XIS external triggerX-ray filename Uniform-Vertical2 ☐ XIS - file closed

Lateral X-ray #2

Time (HH:MM) 11:55

breathhold: _____ secs

Exposure: 77 kV 50 mA 1000 ms☒ Detector position ☒ Trigger cable set ☒ XIS external triggerX-ray filename Uniform-lateral2 ☐ XIS - file closedVector: X 0 Y 0 Z -1 Position: X 58 Y 31.4 Z -11.4

Vertical X-ray #3

Time (HH:MM) 12:16

breathhold: _____ secs

Exposure: 77 kV 50 mA 1000 ms☒ Detector position ☒ Trigger cable set ☒ XIS external triggerX-ray filename Uniform-Vertical3 ☐ XIS - file closed

Lateral X-ray #3

Time (HH:MM) 12:13

breathhold: _____ secs

Exposure: 77 kV 50 mA 1000 ms☒ Detector position ☒ Trigger cable set ☒ XIS external triggerX-ray filename Uniform-lateral3 ☐ XIS - file closed

Vector: X _____ Y _____ Z _____ Position: X _____ Y _____ Z _____

Vertical X-ray #4

Time (HH:MM) _____

breathhold: _____ secs

Exposure: _____ kV _____ mA _____ ms

☐ Detector position ☐ Trigger cable set ☐ XIS external triggerX-ray filename _____ ☐ XIS - file closed

Lateral X-ray #4

Time (HH:MM) _____

breathhold: _____ secs

Exposure: _____ kV _____ mA _____ ms

☐ Detector position ☐ Trigger cable set ☐ XIS external triggerX-ray filename _____ ☐ XIS - file closed

Vector: X _____ Y _____ Z _____ Position: X _____ Y _____ Z _____

Pig Uniform

Date: 24. Juli 2014

Start Field 1: _____ End Field 1: _____

Respirator Log: respirator-log-2014-07-24-12-28-28.txt

Beckhoff Log: Uniform-90.srd

Bemerkungen: sehr lange Breathholds
ein paar "Mini-atmer", PET-Messung nach
der Bestrahlung

Start Field 2: _____ End Field 2: _____

Respirator Log: respirator-log-2014-07-24-13-25-23.txt

Beckhoff Log: Uniform-270.srd

Bemerkungen: anders, als sonst ein paar Mal
während des Spillgate zwischen Respirate ON/OFF
geschaltet

Checkliste Bestrahlungsplanung

DOSE: 40Gy TARGET: PV SCHWEIN: Uniform

Scripts path: AIXd/user/motion/Beamtime/GSI1407/Simulations/SCRIPTS

Heartbeats during CT: N:45 bpm Finalmotion: sin_3mm_1200ms0

| Task | Script | Done? | Name | Date |
|---------------------------------------|--------|-------|------|--------|
| Create Header | | ✓ | AE | 17.07. |
| Copy CT data to PatientData2 | | ✓ | AE | 17.07. |
| Sort CTs (contrast/native) | | ✓ | AE | " |
| Sort motion phased | | ✓ | AE | " |
| Check # of files in each phase | CTX.sh | ✓ | AE | " |
| Check contrast and native slices | " | ✓ | AE | " |
| Create CTX | " | ✓ | AE | " |
| Create MHA | " | ✓ | AE | " |
| Check BB coordinates and write Header | | ✓ | AE | " |
| Copy Contours | | ✓ | AE | 18.07. |
| Create VDX | | ✓ | AE | " |
| Change 0 to 1 for Target (VDX) | | ✓ | AE | " |
| Check voi names in VDX | | ✓ | AE | " |

| | | | | |
|--|--|---|----|--------|
| FalseColorCheck contrast<->native (00) | | ✓ | AE | 17.07. |
| Contours native ok? | | | | |
| Registration contrast -> 4DCT | | ✓ | AE | 17.07. |
| Registration ok? | | ✓ | AE | 17.07. |

| | | | | |
|---------------------------|-------------|---|----|--------|
| CTV - add Margins | Pig_DATA.sh | ✓ | AE | 18.07. |
| CTV - calculate Isocenter | " | ✓ | " | 18.07. |
| ITV - created | " | ✓ | " | 18.07. |
| ITV - calculate Isocenter | " | ✓ | " | 18.07. |

| | | | | |
|---------------------|---|--|--|--|
| For Abl_LV Target: | | | | |
| Create Skin Kontour | " | | | |

| | | | | |
|------------------------------------|-----------------|---|----|--------|
| Static treatment planning | Pig_PLAN.sh | | | |
| Water plan | " | | | |
| Rescanning treatment plannint | " | | | |
| GD Plots | Pig_DOSECUTS.sh | | | |
| DVHs Target + OAR | Pig_ANALYSE.sh | | | |
| Max-Mean Dose Plots | " | | | |
| Dose Volume Plots | " | | | |
| Schweinchenzettel | Pig_FINAL.sh | | | |
| Copy of dose+rst+gd in results dir | " | | | |
| Create MIP and MEAN | Pig_MIP_MEAN.sh | ✓ | AC | 19.07. |

| | | | | |
|----------------------------------|--|--|--|--|
| angle changed from 270° to -90°? | | | | |
| Check BB Isocenter in DY | | | | |
| OVERALL APPROVAL | | | | |

Study Documentation Pig

Uniform

Study group: ☐ AV ☒ PV ☐ LVDose 40 Gy**Surgery:** Date (DD.MM.YYYY): _____

Duration _____ h

Procedures ☐ Pacemaker ☐ ClipsMapping: ☐ LV ☐ RA ☐ LA

Comments: _____

Follow-up #1: Date (DD.MM.YYYY): _____

Duration _____ h

Comments: _____

Follow-up #2: Date (DD.MM.YYYY): _____

Duration _____ h

Comments: _____

CT:Date (DD.MM.YYYY): 16.07.14

Duration _____ h

Scans: ☒ BPL-CT ☐ 4D-CT ☐ 4D-CT with CAHIT patient name: Schwein - 2014 07 16, UniformComments: 38.5 kg**Treatment planning:**

RST filename: _____

Comments: _____

Verification: Date: _____ Time (HH:MM) _____

Dose files: 3D: _____ 4D: _____

Comments: _____

Irradiation: Date: _____ Time (HH:MM) _____ Duration _____ h

SAM filename: _____

Lmdout filename: _____

Comments: _____

CT Documentation Pig Uniform**Premedication:** Ketamin [mg] 150 Azaperon [mg] 200**Medication:** Propofol [ml/KG/min] 0.18 Dormicum [mg] 15Propofol [ml] 180.6

Other: _____

CT: Date (DD.MM.YYYY): 16.07.14 Duration _____ h 18:15 - 19:26HIT patient name: Schwein - 20140716, Uniform ☒ table height 230**Immobilization:** (use room laser on mask) PiggyFix [A/B] AMask Index left: 4-8 right: 4-8Edding: ☒ Laser ☒ ECG lines ☒ Beekleys on laser intersections ☒ Pig nameTattoos: ☒ Laser / mask intersection ☒ mask borders ☒ edge of piggyfixPictures: ☒ left ☒ right ☐ transfer for positioningComments: _____
_____Table position zeroed with CT-laser on 1st frame ☒ other: _____**Respiration:** ☒ 14 bpm, IPPV, V_t = 350, PEEP 0☐ other: _____**Breath-hold test** before 4DCTs ☒ Native ☐ Contrast**Native 4D-CT** Time (HH:MM) 19:06 breathhold: 35 secs# slices: 1331 CTDI_{vol} 63.39 DLP 820.3 Heart rate [bpm]: 46**Contrast agent:** ☒ air-free Dose: ☒ 50 ml ☐ other: _____**Contrast 4D-CT:** Time (HH:MM) 19:08 breathhold: 42 secs# slices: 1331 CTDI_{vol} 63.39 DLP 820.0 Heart rate [bpm]: 50Delay after CA: ☒ 10 sec ☐ other: _____Deviations from protocol "Schwein_GSI": ☐ No ☐ Yes (details in comments!)Comments: 195 mA, (native + contrast)

_____**4D-CT Check** by _____ ☐ Native ☐ Contrast**Contours** by _____ Date (DD.MM.YYYY): _____**Data Transfer to GSI:**☐ 4D-CT ☐ CA-4D-CT Date (DD.MM.YYYY): _____☐ Contours Date (DD.MM.YYYY): _____